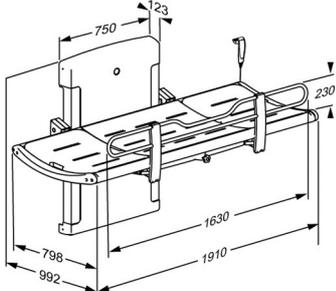




<p><b>Product</b> <b>Item number</b></p>	<p>SCT 3000 shower change table, electrically height adjustable R8582318021</p>
<p><b>Colour</b>  <b>Product description</b></p>	<p>VVS: 777543494 EAN: 5708590357515 HMI: 125822</p> <p>Lead grey</p>  <p>SCT 3000 shower change table, powered, height adjustable 700 mm; from 300 mm - 1000 mm, with wired hand control Foldable. Incl. safety rail and water collection tray with integrated outlet. Surface length 1800 mm, outer measurements 1057 mm x 1910 mm. Excl. mounting kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types. <i>Replaces R8518 &amp; R8538.</i></p> <p>SCT 3000 shower change table, powered, height adjustable 700 mm; from 300 mm - 1000 mm, with wired hand control Foldable. Incl. safety rail and water collection tray with integrated outlet. Surface length 1800 mm, outer measurements 1057 mm x 1910 mm. Materials: Powder-lacquered stainless steel, aluminium, polyurethane foam and polystyrene. Excl. mounting kit.</p>
<p><b>Adjustability</b></p>	<p>Height adjustable</p>
<p><b>Adjustment height</b></p>	<p>700 mm</p>
<p><b>Product capacity</b></p>	<p>Maximum load: 200 kg. When the table is loaded with 200 kg, each mounting bolt supports 26 kg <i>This product has been tested statically with a load of 500 kg as the IEC 60601-standard requires testing with safety factor 2.5.</i></p>
<p><b>Speed</b></p>	<p>Speed up mm: 15 mm/s Speed down mm: 20 mm/s</p>
<p><b>Noise level</b></p>	<p>48 dB</p>
<p><b>Materials</b></p>	<p>Aluminium, polyurethane foam, plastics, powder-lacquered steel.</p>





<b>Product</b>	SCT 3000 shower change table, electrically height adjustable
<b>Item number</b>	R8582318021
<b>Ergonomy</b>	<p>Pressalit offers a range of high quality, fixed or power adjusted shower and changing tables designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centres, institutions, hospitals and private homes.</p> <p>The shower and changing tables 3000 series provides for maximum safety and comfort. Besides an award winning design, the series includes comfortable and durable surfaces, including polyurethane foam and aluminum.</p> <p>The flexibility of the shower and changing table provides for optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p>
<b>Features</b>	<p>A water collection tray with integrated water outlet prevents water from splashing on carers or the floor. The flexible hose attached to the centre mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall. Adjustable neck support which can be moved to either end of the 1800 mm shower and changing tables.</p> <p>The slightly curved form makes the table extremely comfortable and safe.</p> <p>The safety rail can be folded easily up or down.</p> <p>The integrated gas cylinders enable the table to be folded out of the way when not in use to maximise floor space.</p>
<b>Mains connection (V)</b>	230
<b>Actuator (motor)</b>	LA40CS Careline (Linak) 4000 N, IPX6
<b>Control box</b>	SMPS14 Careline (Linak), IPX6
<b>Wired remote control</b>	HB3 (Linak), IPX6
<b>Power consumption (maximum load)</b>	Power consumption by maximum load: 144 W Power consumption by standby: 2 W
<b>Warranty</b>	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is installed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product. The damaging effects can, however, be reduced if the product is showered with tap water on a daily basis.</p>





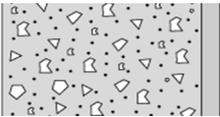
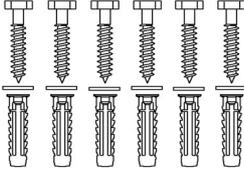
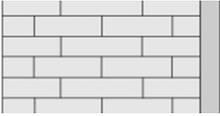
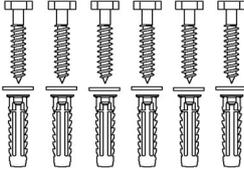
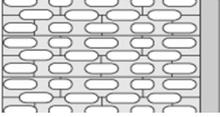
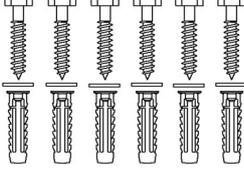
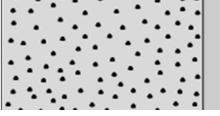
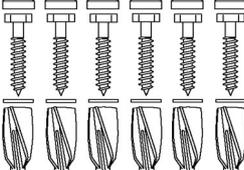
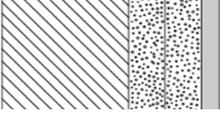
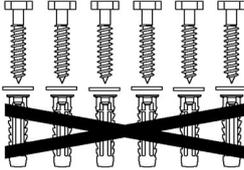
<b>Product</b>	SCT 3000 shower change table, electrically height adjustable
<b>Item number</b>	R8582318021
<b>Cleaning</b>	<p>Showering and changing table 3000. The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Clinell Universal Wipes Clinel Sporidical Wipes Chlor-Clean Tablets WipeClean Ethanol Disinfection 80% Wetwipe Ethanol Disinfection 70% (wetwipe)</p> <p><b>NON-DISINFECTANTS:</b> Wetwipe Universal</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
<b>WEEE</b>	Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.
<b>Declaration of conformity</b>	This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007. It also complies with all relevant requirements in the Council Directive 2006/42/EC (of 17 May 2006) on machinery and the Council Directive 2011/65/EC on RoHS.





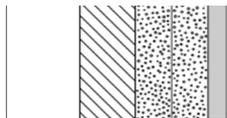
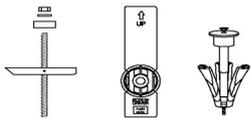
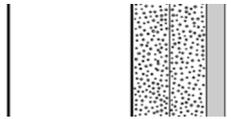
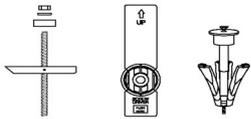
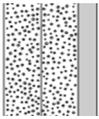
<b>Product</b>	SCT 3000 shower change table, electrically height adjustable
<b>Item number</b>	R8582318021

## Installation in different wall types

	<p><b>Concrete &gt;B25</b></p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p><b>Solid Brick &gt;= Mz12</b></p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p><b>Solid Sand Lime Brick &gt;= Hlz12</b></p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p><b>Aerated Concrete P4</b></p> <p>V8658 - Mounting kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.</p>	
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b></p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.</p>	





<b>Product</b> <b>Item number</b>	SCT 3000 shower change table, electrically height adjustable R8582318021	
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b> Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm)</b> Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	
	<b>Drywall (2 x 12,5 mm) water-resistant</b> Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.	

It is the responsibility of the installer to ensure that suitable wall fasteners are used.



[www.pressalit.com](http://www.pressalit.com)